
BLS CPR Pro SM201

Highest Quality BLS CPR Training

Proven to Achieve Better Quality CPR

A high fidelity manikin with a life-like look and feel. It surpasses other systems in accuracy of measurement, feedback and reporting on more parameters that matter to patients.

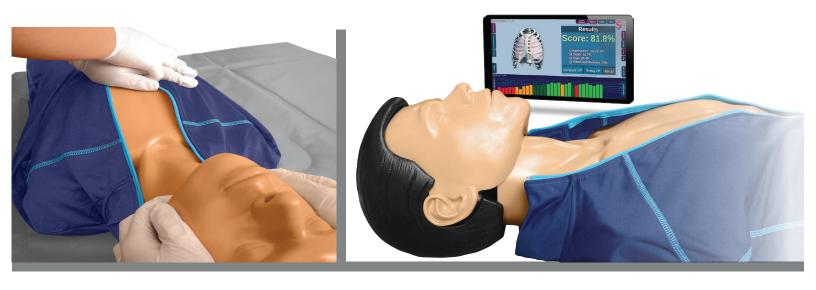

The Only System Giving Accurate Complete Feedback

This system offers real depth and simple to adjust chest resistance for more realistic practice. It detects leaning on the chest and is sensitive to movement away from the correct spot. Ventilations provide accurate display of volume and accurate real time display for rate as you provide air.

Jaw Thrust

iii ust

Metronome

High Quality Feel on Compressing

External Change for Chest Resistance

See As Air Flows Into the Lungs

Natural Resistance on Ventilating



With Cardio-Cerebral Resuscitation

Cerebral Perfusion

Responds to quality and timing of skills performance. Updated each second.

Tidal Flow

Improves use of the BVM. You see accurate flow of air with Volume and Rate as it moves into the lungs.

16155 Sierra Lakes Parkway Suite 160210, Fontana, CA USA 92336 Tel: +1 877 774-4232



Note: Sold as Torso Manikin. Arms & Legs, and IV Training Arm available.

Features Sheet

BLS CPR Pro SM201

High Quality Features Accurate and Complete

AWARD WINNING

Real Time Feedback

- At-a-Glance Understanding of Colored Bars
- More Accurate Than the Competition
- Covers More Parameters

BLS CPR Skills (AHA & ERC Guidelines 2020)

- Compressions: 30:2 x 2 cycles
- Compressions: 30:2 x 4 cycles
- Ventilations: Rescue Breaths
 (feedback on volume, rate, & Interval)
- CPR: 1 AND 2 Rescuer (practice and test modes)

HIGHEST QUALITY

Quick & Easy To Use

- Immediate percentage score
- Display features at a single click
- Display results immediately or later
- Display chart of results
- Single click to print results

Interface

- Compression accuracy displayed in Real Time (Depth, Rate, Recoil, Too Deep, Hitting Bottom)
- Ventilation accuracy displayed in Real Time
- CPR accuracy displayed in Real Time
- Turn Feedback On or Off
- With Feedback off Results Display When Finished

OUALITY CONTROL

Quality Control

- Create User Name (password protected)
- Saves Comments with Results
- Save and Print Results with Names
- Print Class Lists

Advanced Training with Simulations and Trainer Tools

MORE FEATURES

Metronome

- Turn On / Off Metronome
- Defaults to Rate of 110 per minute
- User Selectable Rate

Use With External Screen

- SVGA or HDMI
- Use with External Monitor or Projector

Customize

- Set Number of Cycles
- Set Number of Compressions in Each Cycle

Quick Print Class List

Assemble Lists of Results –
 with Names, Scores, Dates (On Local PC)

SPECIAL ACTIVITIES

Ventilations Trainer FOR 30:2

- Displays Volume and Rate
- Displays Interval & Total Time for the Two Vents

BVM Trainer

- See Rate of Air Flow Into Lungs
- See if BVM is Correctly Released
- See if Volume is Correct

EXTRAS

Ventilation Feedback

- Visually See Volume as it is Delivered
- Improve Technique for Squeezing BVM
- Improve Volume & Rate Control
- Learn Correct Tidal Control of Breathing

Cerebral Perfusion

- Responds in Real Time as Skills are Performed
- Displays a Color Graph of Fraction
- Shows Influence on Blood Flow

